



Littmann®
Stethoscopes

3M™ Littmann® Stethoscopes let you **listen like no other.**

Hear subtle sounds more than twice as loud as the next leading stethoscope.*

There's a reason doctors, nurses, clinicians and healthcare students just embarking on their careers know that 3M™ Littmann® Stethoscopes are synonymous with professionalism, excellent acoustics and dependability. It's because quality is at the core of what we do, and when it comes to choosing a stethoscope, there's no comparison. There's a name in auscultation you can count on—**just look for the L.**

Hear the Hard Stuff

3M™ Littmann® Cardiology IV™ Stethoscopes and 3M™ Littmann® Classic III™ Stethoscopes are more than twice as loud as the next leading stethoscope, allowing you to more reliably and consistently hear heart sounds at lower frequencies (below 120Hz) like Korotkoff sounds, Mitral Stenosis and S3 and S4 Gallops when compared to other leading stethoscopes.*

Handcrafted Quality

Littmann Stethoscopes are designed and manufactured in the United States with globally sourced materials following rigorous processes and strict controls. The result? Littmann stethoscopes are reliably louder than other stethoscopes.*

Confidence and Comfort

Littmann stethoscopes let you confidently assess, monitor and diagnose patients with outstanding acoustics that help you hear subtle sounds other stethoscopes may miss. They are also lighter than other stethoscopes**, providing enhanced all-day comfort.

Most recent innovation:

3M™ Littmann® CORE Digital Stethoscope

- Up to 40x amplification (at peak frequency vs. analog mode)
- Active noise cancellation
- In-app sound wave visualization when used with the Eko app
- Toggle between analog and amplified listening modes to suit your needs
- Plus other powerful features designed to help give you a clearer picture of the health of your patients

Learn more at Littmann.com/CORE




The difference is in the details.

Every 3M™ Littmann® Stethoscope is precision manufactured with advanced materials and features a high level of quality that sets them apart. From top to bottom, here is what makes Littmann stethoscopes stand out.



Lightweight Design

Littmann stethoscopes are designed to feel light and comfortable when draped around your neck or carried in a pocket, giving you an extra step in your day.

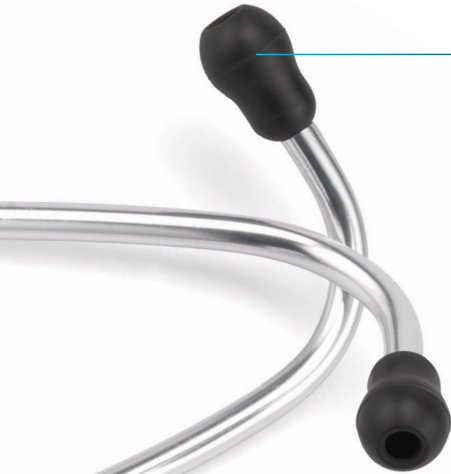

Snap Tight Soft-Sealing Eartips

The bulbous ends of the eartips conform to individual ears for an excellent acoustic seal and comfortable fit. Be sure to use the right size eartips for the best seal.

Adjustable Headset

Headset tension can be easily adjusted for individual fit and comfort. Angled eartubes align with ear canals for comfort and maximum sound transmission.

High Quality, Durable Tubing

Kink-resistant flexible tubing that is easy to clean with 70% isopropyl alcohol—helping to keep your Littmann stethoscope looking new for years to come.

Tunable Diaphragm Technology

Hear different frequency sounds without turning over and repositioning the chestpiece. Just hold the chestpiece with light pressure to hear low frequency sounds, or press a little more to hear higher frequency sounds.


Precision Chestpieces

Designed to provide excellent sound quality over a lifetime of use, chestpieces are crafted with stainless steel and high-density zinc alloy due to their excellent acoustic properties and their strength to resist impacts, scratches and chemicals. Double-sided chestpieces have precision stems that rotate smoothly and click into position.



3M™ Littmann® Cardiology IV™ Stethoscopes and 3M™ Littmann® Classic III™ Stethoscopes

Globally-available comparable stethoscopes in an equivalent class

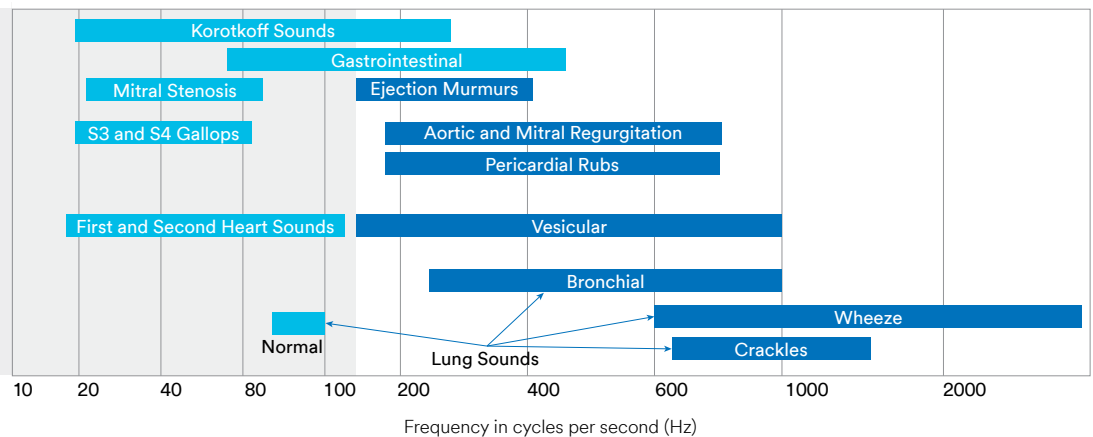

Auscultatory Sounds

Compared to other leading stethoscopes, Littmann stethoscopes excel at allowing you to hear sounds at lower frequencies.*

To learn more, visit www.Littmann.com



3M Medical Solutions Division
2510 Conway Ave.
St. Paul, MN 55144 USA

Listen first.

1-800-228-3957 | www.littmann.com

*Based on tests of the 3M™ Littmann® Cardiology IV™ Stethoscope and 3M™ Littmann® Classic III™ Stethoscope against globally-available comparable stethoscopes in an equivalent class using recorded heart sounds with diaphragm.

**Based on published weights of globally-available comparable stethoscopes in an equivalent class.

3M, Littmann, the L Littmann logo, the L, Cardiology IV and Classic III are marks and/or registered marks of 3M. Unauthorized use prohibited. Used under license in Canada.

Eko trademarks, including CORE, are marks and/or registered marks of Eko Devices Inc.

Please recycle. Printed in USA. © 3M 2024. All rights reserved. 70-2013-1681-0